



NHP  
RESEARCH  
ALLIANCE

---

# **bKIT *Euterpe oleracea***

Real-Time PCR assay

Code: bKTB-EO.01

---

## **Hyris Ltd**

---

### **Hyris Headquarters**

Lower Ground Floor, One George Yard,  
EC3V 9DF, London UK  
Phone: +44.2036082968  
Mail: office@hyris.net

---

### **Hyris Research Center**

Corso Garibaldi, 60  
20121 Milano, Italy  
Phone: +39.02.82951302  
Mail: administrator@hyris.net

---

### **Hyris Asia Pac**

38 Ang Mo Kio Industrial Park 2 #02-07A  
569511 Singapore  
Phone: +65.8160.7207  
Mail: office@hyris.net

---

[www.hyris.net](http://www.hyris.net)

## Hyris Ltd

### Hyris Headquarters

Lower Ground Floor, One George Yard,  
 EC3V 9DF, London UK  
 Phone: +44.2036082968  
 Mail: office@hyris.net

### Hyris Research Center

Corso Garibaldi, 60  
 20121 Milano, Italy  
 Phone: +39.02.82951302  
 Mail: administrator@hyris.net

### Hyris Asia Pac

38 Ang Mo Kio Industrial Park 2 #02-07A  
 569511 Singapore  
 Phone: +65.8160.7207  
 Mail: office@hyris.net

www.hyris.net

#### WARNING AND PRECAUTIONS

Please check Kits integrity before use. Use of deteriorated Kits may cause lack of results and/or equivocal results. Do not mix the components of kits with different batch number. Do not interchange reagents with other kits, or components from different lots. Do not use reagents after their expiration date. This kit is designed to be used by personnel trained to follow correct molecular biology precautions.

#### WARRANTY AND RESPONSIBILITIES

Hyris Ltd guarantees the buyer exclusively the quality of reagents and of the components used to produce the Kits. Any product not fulfilling the specifications included in the product sheet will be replaced. This warranty limits our responsibility to the replacement of the product. No other warranties, of any kind, express or implied, including, without limitation, implicit warranties of commercialisation ability or adequacy for a given purpose, are provided by Hyris Ltd. Hyris Ltd is not responsible and cannot anyway be considered responsible or jointly responsible for possible direct and indirect damages resulting from the utilization of the Kits by the user. The user consciously and under his/her own responsibilities decides for the utilization purposes of the Kits and uses it the way he/she considers most suitable in order to reach his/her goals and/or objectives. Hyris Ltd will not be held responsible for any direct, indirect, consequential or incidental damage resulting of the use, misuses, results of the use or inability to use any product. Some of the applications which may be performed with this product may be covered by applicable patents in certain countries. The purchase of this product does not include or provide a license to perform patented applications. Users may be required to obtain a license depending on the country and/or application. Hyris Ltd does not encourage the unlicensed use of patented applications. The Kits may require the use of Taq Polymerase enzyme, DNA binding components and fluorochromes/quencher, often registered as trademark by companies. TaqMan® is a trademark of Roche Molecular Systems, Inc. FAM™, VIC®, HEX™ and ROX™ are a trademark of Thermo Fisher Scientific or its subsidiaries. SYBR is a registered trademark of Molecular Probes, Inc. (Thermo Fisher Scientific). MGB probes are a trademark of Elitech group. The Kits have been internally tested by our quality control. Any responsibility is waived if the warranty of quality control does not refer to the specific Kits. The user is personally responsible for data that he/she will obtain and/or he/she will supply to third parties using these kits. Once the sealed package is opened the user accepts all the conditions without fail; if the package is still sealed the kit can be returned and the user can be refunded. Kits components are intended, developed, designed, and sold for Research Purpose Only. Product claims are subject to change. Therefore, please refer to our website (www.hyris.net) for the most up-to-date information on Hyris Ltd products.

## *Euterpe oleracea*

*Euterpe oleracea*, commonly known as acai, is a plant species belonging to Arecaceae family. It is widely diffused in central and south America where it has been used since ancient time as food and for treatments of various ailments. Modern analyses confirmed the presence of bioactive compounds, such as anthocyanins and other antioxidant compounds<sup>(1)</sup>.

<sup>(1)</sup> 1: Ulbricht C, Brigham A, Burke D, Costa D, Giese N, Iovin R, Grimes Serrano JM, Tanguay-Colucci S, Weissner W, Windsor R. An evidence-based systematic review of acai (*Euterpe oleracea*) by the Natural Standard Research Collaboration. *J DietSuppl.* 2012 Jun;9(2):128-47. doi: 10.3109/19390211.2012.686347. Review. PubMedPMID: 22607647.

## Principle

Hydrolysis probe Real-Time PCR (qPCR) assay for the detection of *Euterpe oleracea*. The product is intended for research purpose only.

## NHPRA validation

In the validation trials performed by NHPRA (Natural Health Product Research Alliance) the following species were tested: *Serenoa repens*, *Olea europaea*, *Actaea racemosa*, *Beta vulgaris*, *Manihot esculenta*, *Triticum aestivum*, *Vitis vinifera* and *Zea mays*.

## bKIT *Euterpe oleracea* packaging

### Part Number: bKTB-EO.01 -50

|                                      |          |
|--------------------------------------|----------|
| qPCR Master Mix (1 tube, blue cap)   | 50 tests |
| Positive Control (1 tube, green cap) | 10 tests |
| Negative Control (1 tube, red cap)   | 10 tests |

### Part Number: bKTB-EO.01 -100

|                                      |              |
|--------------------------------------|--------------|
| qPCR Master Mix (2 tubes, blue cap)  | 2 x 50 tests |
| Positive Control (1 tube, green cap) | 20 tests     |
| Negative Control (1 tube, red cap)   | 20 tests     |

## Storage

-20°C. Avoid prolonged exposure to light and repeated freeze and thaw cycles.

## Shelf life

If the bKIT is correctly stored, at constant-temperature freezer, its performance is guaranteed until the shelf life indicated on the tubes.

## Additional material/reagents required

- DNA extraction tools and reagents.
- Nuclease-free water.
- Gloves.
- Pipettes.
- bCUBE® instrument or other Real-Time PCR instrument (\*) with filters calibrated for FAM.
- bCUBE® sample loading cartridge or, if using other Real-Time PCR instrument, samples loading support according to the instrument specifications.

(\*) *This assay was especially developed to be used in association with the bCUBE® instrument, available from Hyris Ltd, but can be used also with any other compatible thermal cyclers.*

## Hyris Ltd

### Hyris Headquarters

Lower Ground Floor, One George Yard,  
 EC3V 9DF, London UK  
 Phone: +44.2036082968  
 Mail: office@hyris.net

### Hyris Research Center

Corso Garibaldi, 60  
 20121 Milano, Italy  
 Phone: +39.02.82951302  
 Mail: administrator@hyris.net

### Hyris Asia Pac

38 Ang Mo Kio Industrial Park 2 #02-07A  
 569511 Singapore  
 Phone: +65.8160.7207  
 Mail: office@hyris.net

www.hyris.net

## WARNING AND PRECAUTIONS

Please check Kits integrity before use. Use of deteriorated Kits may cause lack of results and/or equivocal results. Do not mix the components of kits with different batch number. Do not interchange reagents with other kits, or components from different lots. Do not use reagents after their expiration date. This kit is designed to be used by personnel trained to follow correct molecular biology precautions.

## WARRANTY AND RESPONSIBILITIES

Hyris Ltd guarantees the buyer exclusively concerning the quality of reagents and of the components used to produce the Kits. Any product not fulfilling the specifications included in the product sheet will be replaced. This warranty limits our responsibility to the replacement of the product. No other warranties, of any kind, express or implied, including, without limitation, implicit warranties of commercialisation ability or adequacy for a given purpose, are provided by Hyris Ltd. Hyris Ltd is not responsible and cannot anyway be considered responsible or jointly responsible for possible direct and indirect damages resulting from the utilization of the Kits by the user. The user consciously and under his/her own responsibilities decides for the utilization purposes of the Kits and uses it the way he/she considers most suitable in order to reach his/her goals and/or objectives. Hyris Ltd will not be held responsible for any direct, indirect, consequential or incidental damage resulting of the use, misuses, results of the use or inability to use any product. Some of the applications which may be performed with this product may be covered by applicable patents in certain countries. The purchase of this product does not include or provide a license to perform patented applications. Users may be required to obtain a license depending on the country and/or application. Hyris Ltd does not encourage the unlicensed use of patented applications.

The Kits may require the use of Taq Polymerase enzyme, DNA binding components and fluorochromes/quencher, often registered as trademark by companies. TaqMan® is a trademark of Roche Molecular Systems, Inc. FAM™, VIC®, HEX™ and ROX™ are a trademark of Thermo Fisher Scientific or its subsidiaries. SYBR is a registered trademark of Molecular Probes, Inc. (Thermo Fisher Scientific). MGB probes are a trademark of Elitech group. The Kits have been internally tested by our quality control. Any responsibility is waived if the warranty of quality control does not refer to the specific Kits. The user is personally responsible for data that he/she will obtain and/or he/she will supply to third parties using these kits. Once the sealed package is opened the user accepts all the conditions without fail; if the package is still sealed the kit can be returned and the user can be refunded. Kits components are intended, developed, designed, and sold for Research Purpose Only. Product claims are subject to change. Therefore, please refer to our website (www.hyris.net) for the most up-to-date information on Hyris Ltd products.

## DNA extraction

Extract DNA from samples according to your usual protocol. If necessary, Hyris Ltd can recommend an extraction method. At this purpose, contact us at [support@hyris.net](mailto:support@hyris.net).

## Reaction set-up

- Thaw all the bKIT components by placing the tubes on ice.
- Gently mix the tubes content by swirling the tubes.
- Spin the tubes to let the content down.
- In new tubes, one for each sample, including the **Negative Control** and the **Positive Control** of the bKIT, prepare the Reaction Mix as shown in the table below:

| Components   | Volume       |
|--|--------------|
| DNA sample or <b>Positive Control</b> or <b>Negative Control</b> | 5 µL         |
| qPCR Mastermix   | 15 µL        |
| <b>Total Volume</b>  | <b>20 µL</b> |

## Cartridge set-up

The procedure described is for the bCUBE® cartridge, but, if using a different Real-Time PCR instrument, the same procedure can be adopted for other loading sample supports with minor modifications.

### 1. Samples set-up

Samples of the following types must be prepared to be loaded on the cartridge:

**Positive Control** for *Euterpe oleracea*.

**Negative Control** for *Euterpe oleracea*.

Sample(s) to be tested.

### 2. Cartridge Loading

- Load the sample prepared as described in the previous section.
- Carefully seal the cartridge with adhesive film in order to avoid any contamination.
- Load the cartridge onto the bCUBE®, then start the run.

## Method set-up

Set up the run method using the following conditions, depending on the instrument you use.

### 1. On the bCUBE®

- Login on the bAPP.
- Set-up “New Analysis” and Select the “*Euterpe oleracea* 1.x” from the “Global recipes” list.
- Specify the “Well types” for each of the loaded sample as follows (**Fig. 1**):  
 “PosCtrl” for the well loaded with *Euterpe oleracea* **Positive Control**.  
 “NegCtrl” for the well loaded with *Euterpe oleracea* **Negative Control**.  
 “Sample” for the wells loaded with samples under analysis.

**Hyris Ltd**

**Hyris Headquarters**

Lower Ground Floor, One George Yard,  
EC3V 9DF, London UK  
Phone: +44.2036082968  
Mail: office@hyris.net

**Hyris Research Center**

Corso Garibaldi, 60  
20121 Milano, Italy  
Phone: +39.02.82951302  
Mail: administrator@hyris.net

**Hyris Asia Pac**

38 Ang Mo Kio Industrial Park 2 #02-07A  
569511 Singapore  
Phone: +65.8160.7207  
Mail: office@hyris.net

www.hyris.net

**WARNING AND PRECAUTIONS**

Please check Kits integrity before use. Use of deteriorated Kits may cause lack of results and/or equivocal results. Do not mix the components of kits with different batch number. Do not interchange reagents with other kits, or components from different lots. Do not use reagents after their expiration date. This kit is designed to be used by personnel trained to follow correct molecular biology precautions.

**WARRANTY AND RESPONSIBILITIES**

Hyris Ltd guarantees the buyer exclusively concerning the quality of reagents and of the components used to produce the Kits. Any product not fulfilling the specifications included in the product sheet will be replaced. This warranty limits our responsibility to the replacement of the product. No other warranties, of any kind, express or implied, including, without limitation, implicit warranties of commercialisation ability or adequacy for a given purpose, are provided by Hyris Ltd. Hyris Ltd is not responsible and cannot anyway be considered responsible or jointly responsible for possible direct and indirect damages resulting from the utilization of the Kits by the user. The user consciously and under his/her own responsibilities decides for the utilization purposes of the Kits and uses it the way he/she considers most suitable in order to reach his/her goals and/or objectives. Hyris Ltd will not be held responsible for any direct, indirect, consequential or incidental damage resulting of the use, misuses, results of the use or inability to use any product. Some of the applications which may be performed with this product may be covered by applicable patents in certain countries. The purchase of this product does not include or provide a license to perform patented applications. Users may be required to obtain a license depending on the country and/or application. Hyris Ltd does not encourage the unlicensed use of patented applications.

The Kits may require the use of Taq Polymerase enzyme, DNA binding components and fluorochromes/quencher, often registered as trademark by companies. TaqMan® is a trademark of Roche Molecular Systems, Inc. FAM™, VIC®, HEX™ and ROX™ are a trademark of Thermo Fisher Scientific or its subsidiaries. SYBR is a registered trademark of Molecular Probes, Inc. (Thermo Fisher Scientific). MGB probes are a trademark of Elitech group. The Kits have been internally tested by our quality control. Any responsibility is waived if the warranty of quality control does not refer to the specific Kits. The user is personally responsible for data that he/she will obtain and/or he/she will supply to third parties using these kits. Once the sealed package is opened the user accepts all the conditions without fail; if the package is still sealed the kit can be returned and the user can be refunded. Kits components are intended, developed, designed, and sold for Research Purpose Only. Product claims are subject to change. Therefore, please refer to our website (www.hyris.net) for the most up-to-date information on Hyris Ltd products.



**Fig 1. Cartridge set-up**

An example of cartridge set-up on the bAPP for one replicate of a sample to be analyzed is shown.


**2. On a compatible Real-Time PCR instrument**

Please, contact us for the protocol set-up on the instrument.

**Reading the results**

**1. On the bCUBE®**

- a. The presence of the target *Euterpe oleracea* in the **Positive Control** or in the **sample** under analysis will generate an amplification curve (**Fig. 2**)

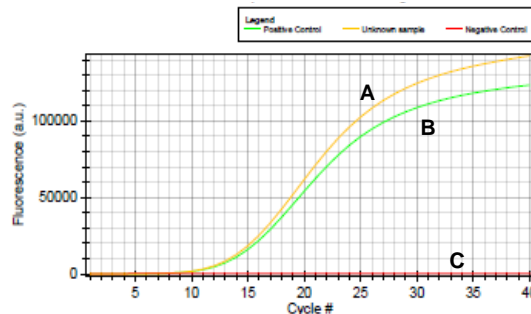

**Fig.2. Amplification plot**

In the plots, the amplification curve of a *Euterpe oleracea* containing **sample (A)**, the **Positive Control (B)**, and the **Negative Control (C)** are shown.

- b. At the end of analysis each well will be labelled depending on the “Well type” as described in the table below and samples classification will be shown on the pdf report of the analysis (**Fig. 3**).

| Well type                  | Possible labels | Label meaning  |
|----------------------------|-----------------|--|
| Positive Control (PosCtrl) | OK              | Specific amplification present                         |
|                            | KO              | Specific amplification absent                          |
| Well type                  | Possible labels | Label meaning  |
| Negative Control (NegCtrl) | OK              | Specific amplification absent                          |
|                            | KO              | Specific amplification present                         |
| Well type                  | Possible labels | Label meaning  |
| Sample                     | Present         | <i>Euterpe oleracea</i> is present in the sample       |
|                            | Absent          | <i>Euterpe oleracea</i> is absent from the sample      |
|                            | Indeterminate   | The test is not conclusive and should be repeated (**) |

(\*\*) If the “Indeterminate” classification persists, contact us at [support@hyris.net](mailto:support@hyris.net).

## Hyris Ltd

### Hyris Headquarters

Lower Ground Floor, One George Yard,  
 EC3V 9DF, London UK  
 Phone: +44.2036082968  
 Mail: office@hyris.net

### Hyris Research Center

Corso Garibaldi, 60  
 20121 Milano, Italy  
 Phone: +39.02.82951302  
 Mail: administrator@hyris.net

### Hyris Asia Pac

38 Ang Mo Kio Industrial Park 2 #02-07A  
 569511 Singapore  
 Phone: +65.8160.7207  
 Mail: office@hyris.net

www.hyris.net

#### WARNING AND PRECAUTIONS

Please check Kits integrity before use. Use of deteriorated Kits may cause lack of results and/or equivocal results. Do not mix the components of kits with different batch number. Do not interchange reagents with other kits, or components from different lots. Do not use reagents after their expiration date. This kit is designed to be used by personnel trained to follow correct molecular biology precautions.

#### WARRANTY AND RESPONSIBILITIES

Hyris Ltd guarantees the buyer exclusively concerning the quality of reagents and of the components used to produce the Kits. Any product not fulfilling the specifications included in the product sheet will be replaced. This warranty limits our responsibility to the replacement of the product. No other warranties, of any kind, express or implied, including, without limitation, implicit warranties of commercialisation ability or adequacy for a given purpose, are provided by Hyris Ltd. Hyris Ltd is not responsible and cannot anyway be considered responsible or jointly responsible for possible direct and indirect damages resulting from the utilization of the Kits by the user. The user consciously and under his/her own responsibilities decides for the utilization purposes of the Kits and uses it the way he/she considers most suitable in order to reach his/her goals and/or objectives. Hyris Ltd will not be held responsible for any direct, indirect, consequential or incidental damage resulting of the use, misuses, results of the use or inability to use any product. Some of the applications which may be performed with this product may be covered by applicable patents in certain countries. The purchase of this product does not include or provide a license to perform patented applications. Users may be required to obtain a license depending on the country and/or application. Hyris Ltd does not encourage the unlicensed use of patented applications.

The Kits may require the use of Taq Polymerase enzyme, DNA binding components and fluorochromes/quencher, often registered as trademark by companies. TaqMan® is a trademark of Roche Molecular Systems, Inc. FAM™, VIC®, HEX™ and ROX™ are a trademark of Thermo Fisher Scientific or its subsidiaries. SYBR is a registered trademark of Molecular Probes, Inc. (Thermo Fisher Scientific). MGB probes are a trademark of Elitech group. The Kits have been internally tested by our quality control. Any responsibility is waived if the warranty of quality control does not refer to the specific Kits. The user is personally responsible for data that he/she will obtain and/or he/she will supply to third parties using these kits. Once the sealed package is opened the user accepts all the conditions without fail; if the package is still sealed the kit can be returned and the user can be refunded. Kits components are intended, developed, designed, and sold for Research Purpose Only. Product claims are subject to change. Therefore, please refer to our website (www.hyris.net) for the most up-to-date information on Hyris Ltd products.

| Results for target <i>Euterpe oleracea</i> |         |
|--|---------|
| Positive control (PosCtrl)                 | OK      |
| Negative control (NegCtrl)                 | OK      |
| Unknown sample (Sample)                    | Present |

Fig.3. Analysis results table

An example of the results table, as reported in the pdf report of the analysis, is shown.

- On a compatible Real-Time PCR instrument  
 Please, contact us for results interpretation.

## Troubleshooting

- Results show no amplification, or anomalous amplification curves

| Possible causes   | Corrective actions   |
|---|--|
| Evaporation of the sample due to inadequate sealing of the plate/strips | Repeat the test using the appropriate materials and tools to seal correctly the plate/strips   |
| Consumables are not appropriate for the method                          | Repeat the test using consumables recommended by the supplier of the Real-Time PCR instrument  |
| The quality of nucleic acid extracted is low                            | Repeat the extraction step. Ensure that the method of extraction has been performed correctly.<br>In any doubt, contact us at <a href="mailto:support@hyris.net">support@hyris.net</a> . |

- No amplification curve is observed for the Positive Control

| Possible causes   | Corrective actions  |
|---|---|
| The Positive Control provided with the assay was not added into the reaction well | Repeat the test adding the Positive Control.<br>If the problem persists, contact us at <a href="mailto:support@hyris.net">support@hyris.net</a> . |

- An amplification curve is observed for the Negative Control

| Possible causes   | Corrective actions   |
|---|--|
| Contamination of the Negative Control or the qPCR Master Mix with target-positive DNA | Repeat the test by applying appropriate quality procedures to prevent contamination.<br>Correctly seal the cartridge or plate/strips.<br>If the problem persists, contact us at <a href="mailto:support@hyris.net">support@hyris.net</a> . |

Document revision Jul 05th, 2019.